

Certol[®]

Medical




*2013 Medical
Product Catalog*



Schedule
Contract GS-07F-236-AA

Developing
the Solutions
For You



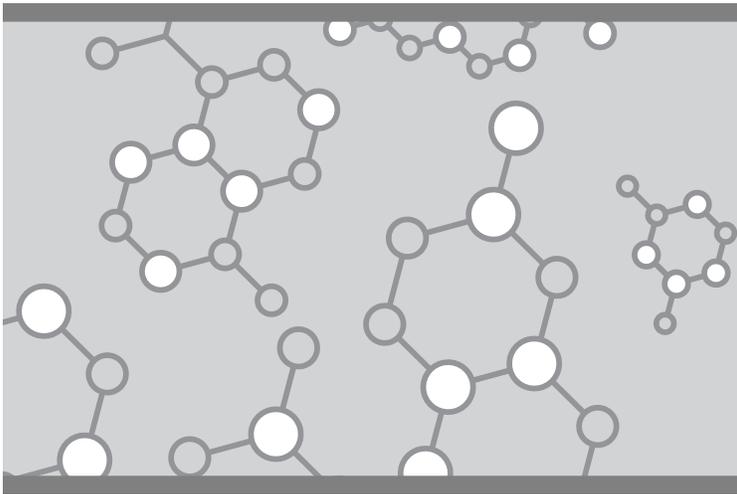

Developing the Solutions For You

Since 1979 Certol® International has been leading the way in the development of infection prevention solutions for you and your facilities. Our internal development and testing ensure that we deliver the very best in cleaning technology for everything from surface disinfection to instrument cleaning and packaging.

Certol is an environmentally responsible manufacturer emphasizing use of biodegradable formulations. Certol follows Good Manufacturing Practices (GMP), and is FDA and EPA compliant. Our production and customer service departments are committed to customer satisfaction with exceptional guarantees.

Our partnerships with organizations like IAHCSSM (International Association of Healthcare Central Service Materiel Management) and SGNA (Society of Gastrointestinal Nurses and Associates) keep us on the cutting edge of what is important to the infection prevention community. We feel that this support of ongoing education is paramount to our continued success.

Certol infection prevention products are in use at major medical centers, ambulatory and surgical centers, doctor's offices, ambulance services, long term care and veterinary clinics. Major medical facilities include the Henry Ford System in Michigan, Cleveland Clinic, Group Health Cooperative in Seattle and many others.

A special thank you to
The University of Colorado Hospital System and
Children's Hospital of Denver
for their cooperation in the production of this catalog.

Surface Disinfectants

Omni II	4
ProSpray C-60	4
ProSpray (Ready-to-Use)	4
ProSpray wipes	5

Enzyme Pre-Treatment

ProEZ foam	6
----------------------	---

Enzyme Instrument Detergents

ProEZ 1	6
ProEZ 2	7
ProEZ AW Quad	7

Instrument/Cart Detergents

ProWash	8
ProClense plus	8

Instrument Lubricants

ProLube Concentrate	9
ProLube RTU	9

Instrument and Cart Rinse Agent

ProRinse	10
--------------------	----

Sterilant Solution

V-Cide	10
------------------	----

Ultrasonic Instrument Detergent

Good Vibrations	11
---------------------------	----

Let the foam Do the Scrubbing!



ProEZ foam™ Ready-to-Use Foaming Enzymatic Spray



Four powerful enzymes (protease, amylase, cellulase and lipase) work to dissolve the broadest range of common surgical soils including proteins, starches, indigestible bioburden and fats.

Certol[®]
Medical

Surface Disinfectants



Omni II™ & ProSpray C-60™ Concentrated Surface Disinfectant/Cleaner

Concentrated disinfectant/cleaner with broad spectrum activity against a wide variety of organisms.

- Effective against MRSA-CA and HA.
- Concentrated format saves money and storage space.
- Approved for 60 day use after dilution.
- Unit dose tubes offer consistent accuracy and convenience. Each tube makes 16 ounces diluted solution.
- Empty 16 oz. spray bottle is labeled for OSHA compliance.



Omni II Description	Case Qty	Order#
32 oz. Tip & Pour Measuring Bottle	12	OMN232-1

ProSpray C-60 Description	Case Qty	Order#
32 oz. Tip & Pour Measuring Bottle	12	PSC60/32-1
Operatory Pack - (48) 1/2 oz. unit dose tubes	3	PSC60/48-1
Intro Kit - (24) 1/2 oz. unit dose tubes, (2) 16 oz. empty spray bottles, (2) squirt caps	1	PSC60/INT-1
16 oz. Empty Spray Bottle	6	PSC60SB-1
Optional Squirt Cap	6	CAPSQUIRT-1

ProSpray™ Ready-to-Use Surface Disinfectant/Cleaner

Disinfectant/cleaner with broad spectrum activity against a wide variety of organisms.

- Effective against MRSA-CA and HA.
- Ready to use format is safe and convenient. No dilution issues!
- Empty 16 ounce bottle is labeled for OSHA compliance. Refill from gallon or 5 gallon drum.
- Water based formula makes it safe for all types of surfaces.
- Fresh lemon scent.



ProSpray Description	Case Qty	Order#
24 oz. Spray Bottle	15	PSC240-1
1 Gallon Refill	4	PSC128-1
5 Gallon Refill with Spigot	1	PSC050-1
Intro Kit - (2) 1 gallon refills, (4) 16 oz. empty spray bottles, (1) funnel, (4) squirt caps	1	PSCSYS-1
16 oz. Empty Spray Bottle	6	PSCPS-1



Kill germs - not equipment with ProSpray disinfectants.

Over time, harsh chemical disinfectants containing alcohols or quaternary ammoniums may degrade upholstery, vinyls and tubing. Peroxides or bleach may cause corrosion of expensive metal fittings and damage color materials.



ProSpray™ wipes Surface Disinfectant/Cleaner Towelettes

- Reliable ProSpray formula with non-alcohol actives. Wipes stay moist longer.
- Protect people AND equipment: Reduce aerosols with ProSpray wipes.
- Ideal for disinfection around electronics, keyboards, digital radiographic sensors and switches.
- Perfect for "high touch" surfaces in public areas implicated in transmission of C. difficile and Norovirus: locker rooms, exercise equipment, TV controls, door knobs, faucets and phones.
- Clean non-submersible power equipment prior to sterilization.



Three Packaging Options:

- Original low-profile soft pack with large 9" x 10" wipe size; stores easily inside drawers.
- Economy canister offers 33 - 50% more wipes than major competitive brands.
- NEW tall high-capacity canister holds 135 extra-large 8-1/2" x 12" wipes. More than DOUBLE the competition.

ProSpray wipes Description	Case Qty	Order#
Soft Pack - 9" x 10" - 72/soft pack	12	PSW-1
Canister - 6" x 6-3/4" - 240/canister	12	PSWC-1
Tall Canister - 8-1/2" x 12" - 135/canister	12	PSWCT-1
Canister Bracket (fits both regular and large canisters)	1	PSWCR-1EA

Water based ProSpray is gentle on surfaces

Great for use on:

- Exam tables
- Wheelchairs
- Positioning cushions
- Stretchers
- Electronics
- OR lights/booms
- Stainless steel tables
- OR stools
- Telephones
- Non-submersibles
- Power drills
- IV poles
- TV remotes
- Bed rails
- Medication control units



ProEZ foam™ Ready-to-Use Foaming Enzymatic Spray

Four powerful enzymes (protease, amylase, cellulase and lipase) work to dissolve the broadest range of common surgical soils including proteins, starches, indigestible bioburden and fats.

- Easy to use – just spray foam over instruments.
- Foam prevents blood soils from drying on instruments while suppressing unpleasant odors.
- Neutral pH formula is gentle on all types of instrumentation.
- Unique anti-corrosive system protects delicate items.
- Eliminate the weight and spills associated with liquid pan presoaks.



See how quickly and effectively ProEZ foam works on tough surgical soils at: www.Certol.com/ProEZfoam

Let the foam do the scrubbing!

ProEZ foam Description	Case Qty	Order#
24 oz. Spray Bottle	15	PREZF240-1
1 Gallon Refill	4	PREZF128-1



ProEZ 1™ Single Enzymatic Detergent

Protease enzymes and detergent are effective on the most common blood based surgical soils.

- Economical yet highly effective.
- Concentrated cleaning power.
- Low foaming formula is versatile and may be used manually, in ultrasonic instrument cleaners and automated washers.
- Neutral pH with anti-corrosive system makes ProEZ 1 safe for use on all types of metals and endoscopes.
- Chelating and rinsing agents improve action in hard water.
- Economical, replacement order option without pump.



Powerful single enzyme detergent attacks the most common protein soils.

ProEZ 1 Description	Case Qty	Order#
1 Gallon with Pump	4	EZ128-1
1 Gallon Re-Order w/o Pump	4	EZNP128-1
5 Gallon Drum	1	EZ050-1
15 Gallon Drum	1	EZ150-1

Enzyme Instrument Detergents



ProEZ 2™ Dual Enzymatic Detergent

Concentrated true dual enzymatic formula includes the concentration and specific enzyme types needed to dissolve complex soils: protease enzymes for blood soils, amylase enzymes for gastrointestinal starch soils and surfactants to break up sticky fats.

- Low foam formula may be used manually, for endoscope flush/soak, in ultrasonic washers and in automated washers.
- Anti-corrosive system and neutral pH protects delicate instruments.
- Odor suppressants and pleasant scent make ProEZ 2 the preferred choice for cleaning scopes and manual processing.
- Chelating and rinsing agents improve action in hard water.
- Economical, replacement order option without pump.



ProEZ AW Quad™ Quadruple Enzyme Automatic Washer Detergent

Certol's premium enzymatic detergent provides the chemistry needed for difficult soils, complex instrumentation and high impingement automated cleaning equipment.

- Four enzymes attack the entire range of surgery soils: blood proteins, gastrointestinal starches, cellulose and fatty orthopedic soils.
- Chelating and rinsing agents improve action in hard water.
- Formula delivers the concentration of active agents needed for fast, consistent results. Starts working quickly – in as little as one minute in warm water.
- Powerful yet safe: biodegradable, neutral pH, instrument brighteners and unique anti-corrosive system. Use on all types of metals, plastics and rubber goods.
- Penetrating agents improve cleaning action around serrations and box locks.
- Economical, replacement order option without pump.



The most effective true dual enzyme detergent on the market.

ProEZ 2 Description	Case Qty	Order#
1 oz. Unit Dose - 24/box	6	PREZU24-1
1 Gallon with Pump	4	PREZ128-1
1 Gallon Re-Order w/o Pump	4	PREZNP128-1
5 Gallon Drum	1	PREZ050-1
15 Gallon Drum	1	PREZ150-1
55 Gallon Drum	1	PREZ550-1

Non-foaming and neutral pH!

ProEZ AW Quad Description	Case Qty	Order#
1 Gallon	4	PREZAW128-1
1 Gallon Re-Order w/o Pump	4	PREZAWN128-1
5 Gallon Drum	1	PREZAW050-1
15 Gallon Drum	1	PREZAW150-1
30 Gallon Drum	1	PREZAW300-1
55 Gallon Drum	1	PREZAW550-1

Instrument/Cart Detergents



ProWash™ Instrument Detergent and Cart Wash

Premium neutral pH multi-purpose detergent combines a powerful cleaning system with effective anti-corrosives.

- Ideal material compatibility to meet manufacturer's instructions. Safe for use on lensed items, aluminum, anodized aluminum, copper, stainless and carbon steel. Solvent-free cleaning agents recommended for items with gaskets, o-rings and other delicate fittings.
- Formulated for use in hard water conditions. High levels of chelating and free-rinsing agents ensure best cleaning results and prevent re-deposition of soils.
- Versatile low foam detergent is effective for use in manual/soak applications, cart washers, ultrasonic instrument washers and high impingement automatic instrument washers.
- Concentrated and cost-efficient: Dilute at ½ to 1 ounce per gallon depending on water quality.



ProClense™ plus Instrument Detergent

ProClense plus is an economical detergent with effective cleaning agents and anti-corrosive agents to help maintain instrument surface quality.

- Moderate alkalinity meets AAMI recommendations for effective cleaning and material compatibility. Safe for use on stainless and carbon steel, aluminum and copper.
- Low foaming formula is effective for use in automatic instrument washers and cart washers.
- Formulated for use in hard water conditions. High levels of chelating and free-rinsing agents ensure best cleaning results and prevent re-deposition of soils.
- Mild clean scent and no added dyes.



Robust formula for wide variety of applications.

Economical and effective cleaning action.

ProWash Description	Case Qty	Order#
1 Gallon	4	PW128-1
5 Gallon Drum	1	PW050-1
15 Gallon Drum	1	PW150-1
30 Gallon Drum	1	PW300-1
55 Gallon Drum	1	PW550-1

ProClense plus Description	Case Qty	Order#
1 Gallon	4	PCNF128-1
5 Gallon Drum	1	PCNF050-1
15 Gallon Drum	1	PCNF150-1
30 Gallon Drum	1	PCNF300-1
55 Gallon Drum	1	PCNF550-1



ProLube™ Concentrate Instrument Lubricant Concentrate

Water soluble instrument lubricant designed for use in the lubrication cycle of automated washer systems.

- Unlike competitors, ProLube Concentrate is a true solution that never separates.
- Unique mineral oil and silicone-free, water soluble, non-toxic food grade lubricant.
- Permeable to steam and ETO sterilants.
- No rinsing necessary after application.
- Penetrating agents increase lubricant action in hinges and box locks.
- Anti-corrosives, rust inhibitors and surface conditioning agents improve instrument function and surface integrity.
- ProLube will never produce a sticky yellow film inside box locks after hot sterilization cycles, unlike traditional mineral oil based "instrument milk".



ProLube™ RTU Instrument Lubricant Ready-to-Use

ProLube RTU is applied manually after cleaning and before packaging and sterilization. Easy to use spray allows spot treatment of box locks, hinges and jointed items.


- Unique mineral oil and silicone-free, water soluble, non-toxic food grade lubricant.
- No rinsing necessary after application.
- Completely permeable to steam and ethylene oxide sterilants.
- A true solution that never separates – no shaking required.
- No oily or sticky build up on instruments.
- Convenient and easy to apply. Eliminates cross-contamination.



According to ANSI/AAMI ST79 "...Water-soluble instrument lubricants, specifically designed for compatibility with sterilization, may be used... Instrument lubricants containing mineral oil or other oil bases should not be used, except to lubricate the internal mechanism of powered instruments..." [2]

ProLube Concentrate

Description	Case Qty	Order#
1 Gallon	4	PLC128-1
5 Gallon Drum	1	PLC050-1
15 Gallon Drum	1	PLC150-1
30 Gallon Drum	1	PLC300-1
55 Gallon Drum	1	PLC550-1

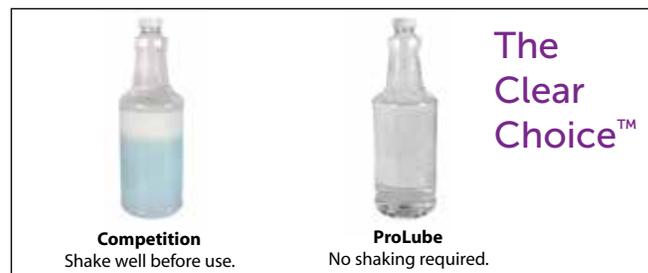

ProLube RTU Description	Case Qty	Order#
32 oz. Spray Bottle	15	PLR/32-1
System Kit - (4) 32 oz. spray bottles, (5) 16 oz. empty spray bottles, (1) funnel	1	PLRSYS-1



ProRinse™ Concentrated Rinse Agent

Concentrated rinse aid for use in automatic instrument and cart washers. Certol's premium clean rinsing agent is designed to reduce drying time and spotting caused by hard water.

- Special formulation designed to reduce drying time.
- Creates sheeting action, leaving carts, glassware and instruments spot free.
- Concentrated formula makes ProRinse as economical as it is effective.



ProRinse Uses:

- Automated washers
- Cart washers
- Mechanical cleaning equipment

Dry faster and spot less with ProRinse.

ProRinse Description	Case Qty	Order#
1 Gallon	4	PRR128-1
5 Gallon Drum	1	PRR050-1
15 Gallon Drum	1	PRR150-1
55 Gallon Drum	1	PRR550-1

V-Cide® Chemical Vapor Sterilant Solution

V-Cide is a chemical sterilant formulated for use in Chemiclave® sterilizers.



Reduced Odor

- Improved formula reduces odor-causing agents.

Gentle and Non-Corrosive

- Does not corrode, rust or dull instruments.

Economical

- Cost-effective alternative for your chemical vapor sterilizer.

Safe

- Compatible with materials recommended by manufacturers for sterilization in chemical vapor sterilizers.



Formulated for use in Chemiclave® sterilizers.

Chemiclave is a registered trademark of Thermo Scientific.

V-Cide Description	Case Qty	Order#
Liter Bottle	16	VC338-1



Good Vibrations™ Ultrasonic Cleaner Solution


Combines great cleaning action with excellent anti-corrosives and instrument brighteners.

- Formula is effective for both organic (blood) soils and inorganic materials.
- Non-foaming formula supports ultrasonic action.
- Safe for delicate surgical items.
- Ideal for jointed pliers and forceps.
- Concentrated economy – 1 ounce makes 1 gallon solution.
- Chelating and rinsing agents improve action in hard water.



Clean and protect against corrosion with Good Vibrations.

Good Vibrations Description	Case Qty	Order#
1 oz. Unit Dose - 24/box	6	GVUS24-1
1 Gallon with Pump	4	GVUS128-1

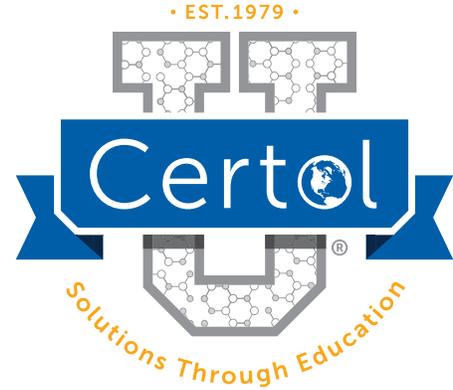

Certol supports quality improvement in SPD with complimentary in-service and continuing education programs to IAHCSSM chapter groups.

Infection prevention is serious business and we understand how everyday activities can challenge even the most dedicated sterile processing professional.

At Certol, we formulate our detergents to work under the toughest of circumstances. Our education programs are designed the same way – to help you navigate difficult cleaning problems and achieve desired “best practices” for your department.

In addition to facility in-service programs for Certol products, Certol Education also offers objective non-commercial programs with experienced speakers on the following topics (speaker CV and program objectives available for review upon request):

- **Manufacturer’s IFU and Material Compatibility:** *Are You Caught in the Middle?*
- **OSHA is Not a 4 Letter Word:** *Chemical Safety in SPD.*
- **Wonderful World of Clean:** *Using Enzyme Science in SPD.*

SCHEDULE TODAY!

Visit: www.Certol.com/University to submit your request.

Programs are offered on a first come, first serve basis, at your facility. Our Clinical Education Department will contact you to schedule a program. Or contact Certol’s Clinical Support Department at (800) 843-3343.

**Offerings based on availability.*

2013 GSA Price List GS-07F-236-AA - Schedule 73



Contract Period: May 1, 2013 to April 30, 2018

Contractor: Certol International, LLC
 6120 E 58th Ave.
 Commerce City, CO 80022-3941
 Ph: 303-799-9401 Fax: 303-799-9408
<http://www.certol.com>

Contractor's Administration Source: Aaron Trapp
atrapp@certol.com

Business Size: Small
Socioeconomic Indicators: None

CUSTOMER INFORMATION:

1a. TABLE OF AWARDED SPECIAL ITEM NUMBERS (SINs)

SIN	DESCRIPTION
375-363	Cleaners/Degreasers
375-372	Cleaning Chemicals Used with Dispensing Systems
476-12	Sterilants, Surgical Scrubs, and Related Products
476-13	Disinfectants

1b. LOWEST PRICED MODEL NUMBER AND PRICE FOR EACH SIN: Prices include freight within CONUS. Volume commitment is required for TIER 2 and TIER 3 pricing.						MSRP	GSA ACCESS PRICE	GSA VOLUME DISCOUNT PRICE (COMMITMENT REQUIRED)	GSA VOLUME DISCOUNT PRICE (COMMITMENT REQUIRED)
FSC Number	Special Item Number (SIN)	Product Part Number	Product Description	UOM	QTY/Case	Suggested Retail Case Price	GSA \$0-4,999 TIER 1	GSA \$5,000-39,999 TIER 2	GSA \$40,000 + TIER 3
79	375-363	PREZF240-1	Pro EZ foam Quadruple Enzymatic Instrument Spray 24 oz. Bottle	CA	15	\$258.75	\$80.50	\$74.74	\$69.76
79	375-363	PREZF128-1	Pro EZ foam Quadruple Enzymatic Instrument Spray Refill	CA	4	\$183.00	\$66.39	\$61.20	\$57.50
79	375-363	EZ128-1	ProEZ 1 Single Enzymatic Detergent Concentrate 1 Gal. w/1 oz. Pump	CA	4	\$128.00	\$59.00	\$54.80	\$51.00
79	375-363	EZNP128-1	ProEZ 1 Single Enzymatic Detergent Concentrate 1 Gal. w/o Pump	CA	4	\$107.00	\$51.31	\$47.64	\$44.47
79	375-372	EZ050-1	ProEZ 1 Single Enzymatic Detergent Concentrate 5 Gallon	DRUM	1	\$140.00	\$73.42	\$68.17	\$63.63
79	375-372	EZ150-1	ProEZ 1 Single Enzymatic Detergent Concentrate 15 Gal.	DRUM	1	\$397.00	\$197.80	\$183.50	\$171.30
79	375-363	PREZU24-1	ProEZ 2 Dual Enzymatic Detergent Concentrate (24) 1 oz. Unit Dose Tubes	CA	6	\$141.00	\$69.35	\$64.40	\$60.11
79	375-363	PREZ128-1	ProEZ 2 Dual Enzymatic Detergent Concentrate 1 Gal. w/ 1 oz. Pump	CA	4	\$141.60	\$64.16	\$59.56	\$55.60
79	375-363	PREZNP128-1	ProEZ 2 Dual Enzymatic Detergent Concentrate 1 Gal. w/o Pump	CA	4	\$122.00	\$54.80	\$50.89	\$47.49
79	375-372	PREZ050-1	ProEZ 2 Dual Enzymatic Detergent Concentrate 5 Gal.	DRUM	1	\$177.00	\$78.00	\$72.43	\$67.60
79	375-372	PREZ150-1	ProEZ 2 Dual Enzymatic Detergent Concentrate 15 Gal.	DRUM	1	\$498.00	\$211.40	\$196.40	\$184.40
79	375-372	PREZ550-1	ProEZ 2 Dual Enzymatic Detergent Concentrate 55 Gal.	DRUM	1	\$1,770.00	\$700.00	\$650.00	\$607.00
79	375-363	PREZAW128-1	ProEZ AW Quad Concentrated Quadruple Enzyme Detergent 1 Gal. w/ 1 oz. Pump	CA	4	\$169.00	\$73.35	\$65.69	\$58.68
79	375-363	PREZAWN128-1	ProEZ AW Quad Concentrated Quadruple Enzyme Detergent 1 Gal. w/o Pump	CA	4	\$137.00	\$55.45	\$51.49	\$48.05
79	375-372	PREZAW050-1	ProEZ AW Quad Concentrated Quadruple Enzyme Detergent 5 Gal.	DRUM	1	\$199.00	\$77.60	\$69.80	\$65.10
79	375-372	PREZAW150-1	ProEZ AW Quad Concentrated Quadruple Enzyme Detergent 15 Gal.	DRUM	1	\$549.00	\$220.33	\$197.31	\$182.20
79	375-372	PREZAW300-1	ProEZ AW Quad Concentrated Quadruple Enzyme Detergent 30 Gal.	DRUM	1	\$1,085.00	\$410.08	\$367.24	\$337.00
79	375-372	PREZAW550-1	ProEZ AW Quad Concentrated Quadruple Enzyme Detergent 55 Gal.	DRUM	1	\$1,968.00	\$740.50	\$663.13	\$626.00

1b. LOWEST PRICED MODEL NUMBER AND PRICE FOR EACH SIN: Prices include freight within CONUS. Volume commitment is required for TIER 2 and TIER 3 pricing.						MSRP	GSA ACCESS PRICE	GSA VOLUME DISCOUNT PRICE (COMMITMENT REQUIRED)	GSA VOLUME DISCOUNT PRICE (COMMITMENT REQUIRED)
FSC Number	Special Item Number (SIN)	Product Part Number	Product Description	UOM	QTY/Case	Suggested Retail Case Price	GSA \$0-4,999	GSA \$5,000-39,999	GSA \$40,000 +
							TIER 1	TIER 2	TIER 3
79	375-363	GVUS24-1	Good Vibrations Multi-Purpose Ultrasonic Cleaner Concentrate - 24 (1 oz.) Unit Dose Tubes	CA	6	\$129.00	\$50.03	\$45.48	\$41.20
79	375-363	GVUS128-1	Good Vibrations Multi-Purpose Ultrasonic Cleaner Concentrate 1 Gal.	CA	4	\$100.00	\$65.62	\$60.00	\$56.00
79	375-363	ICL128-1	ProClense One Step Instrument Detergent/Lubricant 1 Gal.	CA	4	\$70.00	\$31.28	\$28.24	\$25.41
79	375-372	PCL050-1	ProClense One Step Instrument Detergent/Lubricant 5 Gal.	DRUM	1	\$81.75	\$36.51	\$32.96	\$29.66
79	375-372	PCL150-1	ProClense One Step Instrument Detergent/Lubricant 15 Gal.	DRUM	1	\$188.50	\$97.43	\$87.96	\$79.16
79	375-372	PCL300-1	ProClense One Step Instrument Detergent/Lubricant 30 Gal.	DRUM	1	\$325.00	\$167.77	\$151.46	\$136.31
79	375-372	PCL550-1	ProClense One Step Instrument Detergent/Lubricant 55 Gal.	DRUM	1	\$503.50	\$297.12	\$268.24	\$241.41
79	375-363	PCNF128-1	ProClense plus Instrument Detergent 1 Gal.	CA	4	\$70.00	\$31.28	\$28.80	\$26.88
79	375-372	PCNF050-1	ProClense plus Instrument Detergent 5 Gal.	DRUM	1	\$81.00	\$39.32	\$36.51	\$34.08
79	375-372	PCNF150-1	ProClense plus Instrument Detergent 15 Gal.	DRUM	1	\$188.00	\$102.25	\$94.94	\$88.61
79	375-372	PCNF300-1	ProClense plus Instrument Detergent 30 Gal.	DRUM	1	\$325.00	\$174.20	\$161.80	\$151.10
79	375-372	PCNF550-1	ProClense plus Instrument Detergent 55 Gal.	DRUM	1	\$503.00	\$311.00	\$289.00	\$269.60
79	375-363	PW128-1	ProWash Instrument Detergent and Cart Wash Concentrate 1 Gal.	CA	4	\$70.00	\$31.28	\$28.80	\$26.88
79	375-372	PW050-1	ProWash Instrument Detergent and Cart Wash Concentrate 5 Gal.	DRUM	1	\$81.00	\$39.32	\$36.51	\$34.08
79	375-372	PW150-1	ProWash Instrument Detergent and Cart Wash Concentrate 15 Gal.	DRUM	1	\$188.00	\$102.25	\$94.94	\$88.61
79	375-372	PW300-1	ProWash Instrument Detergent and Cart Wash Concentrate 30 Gal.	DRUM	1	\$325.00	\$174.20	\$161.80	\$151.10
79	375-372	PW550-1	ProWash Instrument Detergent and Cart Wash Concentrate 55 Gal.	DRUM	1	\$503.00	\$311.00	\$289.00	\$269.60
79	375-372	PLC128-1	ProLube Concentrate Instrument Lubricant Concentrate 1 Gal.	CA	4	\$97.00	\$61.68	\$57.27	\$53.44
79	375-372	PLC050-1	ProLube Concentrate Instrument Lubricant Concentrate 5 Gal.	Drum	1	\$114.17	\$77.40	\$71.80	\$67.05
79	375-372	PLC150-1	ProLube Concentrate Instrument Lubricant Concentrate 15 Gal.	Drum	1	\$358.00	\$212.25	\$197.09	\$183.95
79	375-372	PLC300-1	ProLube Concentrate Instrument Lubricant Concentrate 30 Gal.	Drum	1	\$742.00	\$407.00	\$378.00	\$353.00
79	375-372	PLC550-1	ProLube Concentrate Instrument Lubricant Concentrate 55 Gal.	Drum	1	\$1,295.00	\$740.00	\$687.00	\$641.00
79	375-372	PLR/32-1	ProLube RTU Ready-To-Use Instrument Lubricant Spray Bottle 32 oz.	CA	15	\$112.50	\$37.99	\$34.53	\$32.15
79	375-372	PLRSYS-1	ProLube RTU Ready-To-Use Instrument Lubricant Sys. Kit: 4 qts, 5 (16 oz.) spray btl, funnel	KIT	1	\$52.00	\$22.18	\$20.58	\$19.22

1b. LOWEST PRICED MODEL NUMBER AND PRICE FOR EACH SIN: Prices include freight within CONUS. Volume commitment is required for TIER 2 and TIER 3 pricing.						MSRP	GSA ACCESS PRICE	GSA VOLUME DISCOUNT PRICE (COMMITMENT REQUIRED)	GSA VOLUME DISCOUNT PRICE (COMMITMENT REQUIRED)
FSC Number	Special Item Number (SIN)	Product Part Number	Product Description	UOM	QTY/Case	Suggested Retail Case Price	GSA \$0-4,999	GSA \$5,000-39,999	GSA \$40,000 +
							TIER 1	TIER 2	TIER 3
79	375-372	PRR128-1	ProRinse Concentrated Instrument and Cart Rinse Agent 1 Gal.	CA	4	\$198.00	\$101.35	\$94.10	\$87.90
79	375-363	PRR050-1	ProRinse Concentrated Instrument and Cart Rinse Agent 5 Gal.	DRUM	1	\$253.00	\$130.00	\$120.70	\$112.60
79	375-363	PRR150-1	ProRinse Concentrated Instrument and Cart Rinse Agent 15 Gal.	DRUM	1	\$723.50	\$372.00	\$345.50	\$322.50
79	375-363	PSW-1	ProSpray wipes Pre-Moistened Surface Disinfectant/Cleaner Soft Pack - 72 (9" x 10") Towelettes	CA	12	\$99.00	\$41.33	\$40.09	\$38.44
68	476-13	PSWC-1	ProSpray wipes Pre-Moistened Surface Disinfectant/Cleaner Canister - 240 (6" x 6-3/4") Towelettes	CA	12	\$117.00	\$47.73	\$46.25	\$44.39
68	476-13	PSWCT-1	ProSpray wipes Pre-Moistened Surface Disinfectant/Cleaner - Canister 135 (8-1/2" x 12") Towelettes	CA	12	\$149.40	\$62.40	\$60.25	\$57.80
68	476-13	PSWCR-1EA	ProSpray wipes Canister Wall Mount Bracket (Branded)	EA	1	\$28.25	\$15.00	\$14.00	\$12.99
68	476-13	PSC240-1	ProSpray Ready-to-use Disinfectant / Cleaner Spray Bottle 24 oz.	CA	15	\$119.25	\$53.63	\$52.02	\$49.88
68	476-13	PSC128-1	ProSpray Ready-to-use Disinfectant / Cleaner Refill 1 Gal.	CA	4	\$96.00	\$27.60	\$26.77	\$25.67
68	476-13	PSC050-1	ProSpray Ready-to-use Disinfectant / Cleaner Refill 5 Gal.	DRUM	1	\$93.00	\$34.50	\$33.45	\$32.05
68	476-13	PSCSYS-1	ProSpray Ready-to-use Disinfectant / Cleaner Kit: 2 (1 Gal.), 4 (16 oz) Empty Spray Btls, Funnel	KIT	1	\$50.00	\$29.45	\$27.35	\$25.55
68	476-13	PSCPS-1	ProSpray Ready to-use-Disinfectant / Cleaner - Empty 16 oz. Empty Pump Spray Bottle	CA	6	\$20.10	\$11.50	\$10.85	\$10.12
68	476-13	PSC60/INT-1	ProSpray C-60 Concentrated Surface Disinfectant/Cleaner - Intro. Kit (24 (1/2 oz.) unit dose, 2 (16 oz.) spray bottles, 2 squirt caps)	KIT	1	\$29.50	\$15.02	\$14.55	\$13.97
68	476-13	PSC60/48-1	ProSpray C-60 Concentrated Surface Disinfectant/Cleaner - Operatory Pack - 48 1/2 oz. unit dose	CA	6	\$330.00	\$54.48	\$52.80	\$50.67
68	476-13	PSC60/32-1	ProSpray C-60 Concentrated Surface Disinfectant/Cleaner 32 ounce Meter Dose Bottle	CA	12	\$478.20	\$117.97	\$114.43	\$109.71
68	476-13	PSC60SB-1	ProSpray C-60 Concentrated Surface Disinfectant/Cleaner- 16 oz. Empty Pump Spray Bottle (with optional squirt Cap)	CA	6	\$20.10	\$11.50	\$11.15	\$10.70
68	476-13	CAPSQUIRT-1	ProSpray C-60 Concentrated Surface Disinfectant/Cleaner - Squirt Cap (Optional)	EA	6	\$0.25	\$0.90	\$0.90	\$0.90
68	476-13	OMN232-1	OMNI II Concentrated Instrument Holding Solution and Disinfectant / Cleaner - 32 ounce Meter Dose Bottle	CA	12	\$351.00	\$160.00	\$155.20	\$148.80
68	476-12	VC338-1	V-Cide Chemical Vapor Sterilant Solution 1 Liter	CA	16	\$148.00	\$70.50	\$65.40	\$61.10

2. MAXIMUM ORDER*: \$150,000 per SINS 375-363, 375-372 and \$125,000 per SIN 476-12, 476-13.

*If the best value selection places your order over the Maximum Order identified in this catalog/pricelist, you have an opportunity to obtain a better schedule contract price. Before placing your order, contact the aforementioned contractor for a better price. The contractor may (1) offer a new price for this requirement (2) offer the lowest price available under this contract or (3) decline the order. A delivery order that exceeds the maximum order may be placed under the schedule contract in accordance with FAR 8.404.

3. MINIMUM ORDER: 1 case

4. GEOGRAPHIC COVERAGE: Domestic, 50 states, Washington, DC, US Territories

5. POINT(S) OF PRODUCTION: Commerce City, CO

6. DISCOUNT FROM LIST PRICES: Basic discount of 31.90%-68.90% from the awarded commercial price list. For calculation of the GSA Schedule price (price paid by customers ordering from the GSA Schedule, and the price to be loaded in to GSA Advantage), deduct the appropriate basic discount from the list price and add the prevailing IFF rate to the negotiated discounted price (Net GSA price).

The current IFF is .75% and should be calculated as follows: Negotiated price divided by (1 minus .0075) which equates to Negotiated price divided by 0.9925. Example: (\$100,000 / 0.9925) = \$100,755.67

7. QUANTITY DISCOUNT(S): 0-13%

8. PROMPT PAYMENT TERMS: Net 30

9.a Government Purchase Cards must be accepted at or below the micro-purchase threshold.

9.b Government Purchase Cards are accepted above the micro-purchase threshold. Contact contractor for limit.

10. FOREIGN ITEMS: None

11a. TIME OF DELIVERY: Shipped 5-7 Days after receipt of order

11b. EXPEDITED DELIVERY: 4-6 days

11c. OVERNIGHT AND 2-DAY DELIVERY: If available, contact the Contractor for rates.

11d. URGENT REQUIREMENTS: Customers are encouraged to contact the contractor for the purpose of requesting accelerated delivery.

12. FOB POINT: Destination CONUS

13a. ORDERING ADDRESS: Same as contractor

13b. ORDERING PROCEDURES: For supplies and services, the ordering procedures, information on Blanket Purchase Agreements (BPA's) are found in FAR 8.405-3

14. PAYMENT ADDRESS: Same as contractor

15. WARRANTY PROVISION: Customer should contact contractor for a copy of the warranty.

16. EXPORT PACKING CHARGES: Not applicable

17. TERMS AND CONDITIONS OF GOVERNMENT PURCHASE CARD ACCEPTANCE: N/A

18. TERMS AND CONDITIONS OF RENTAL, MAINTENANCE AND REPAIR (IF APPLICABLE): N/A

19. TERMS AND CONDITIONS OF INSTALLATION (IF APPLICABLE): N/A

20. TERMS AND CONDITIONS OF REPAIR PARTS INDICATING DATE OF PARTS PRICE LISTS AND ANY DISCOUNTS FROM LIST PRICES (IF AVAILABLE): N/A

20a. TERMS AND CONDITIONS FOR ANY OTHER SERVICES (IF APPLICABLE): N/A

21. LIST OF SERVICE AND DISTRIBUTION POINTS (IF APPLICABLE): N/A

22. LIST OF PARTICIPATING DEALERS (IF APPLICABLE): Cardinal Health; Henry Schein, Inc.; Medline Industries, Inc.; McKesson Medical; Seneca Medical

23. PREVENTIVE MAINTENANCE (IF APPLICABLE): N/A

24a. SPECIAL ATTRIBUTES SUCH AS ENVIRONMENTAL ATTRIBUTES : N/A

24b. Section 508 Compliance for EIT: as applicable

25. DUNS NUMBER: 095153714

26. NOTIFICATION REGARDING REGISTRATION IN SYSTEM FOR AWARD MANAGEMENT (SAM) DATABASE: Contractor has an Active Registration in the SAM database

Customer Satisfaction Guaranteed

Our Customer Service Department will provide you with complete infection prevention methodology and product information and help you achieve perfection in your techniques. Call 1-800-THE-EDGE (1-800-843-3343). Our staff is available weekdays from 7:00am - 4:00pm (MST) to answer your questions.

We can also be reached at www.Certol.com. Our website provides in-depth infection prevention procedures and Material Safety Data Sheets for all Certol products.

Certol is not only committed to providing high quality products and the complete satisfaction of our customers, Certol ABSOLUTELY GUARANTEES IT. If the end user of any Certol product is not completely satisfied with the quality of the product, a refund or adjustment of equal value will be made. For a copy of any of Certol's product warranties call 1-800-THE-EDGE (1-800-843-3343).

Ordering Information

- Product sold in minimum of case quantities only.
- For international orders, please contact Customer Service for specific terms and conditions.

Terms

Net 30 days.

Shipping Information

- F.O.B. CONUS.
- Shipped 5-7 days after receipt of order.
- Contact Customer Service for accelerated delivery.

Return Policy

- A Returned Goods Authorization (RGA) must be obtained from Customer Service prior to a product return.
- Full cases only.
- Customer is responsible for return shipping.
- 25% restocking fee on all returned items.
- 14-month shelf life on products returned.
- 25% disposal fee on unsellable products.
- All product returns must be within six months of receipt of goods. Credit memos issued will expire after six months.

All Certol products are available through national and regional distributors. Please visit: www.Certol.com/MedicalDistributors for the complete list or Call (800) 843-3343 to locate a distributor near you!



Quotations and References:

1. Centers for Disease Control and Prevention, *Guidelines for Disinfection and Sterilization in Healthcare Facilities*, 2008.
2. "Comprehensive guide to steam sterilization and sterility assurance in health care facilities" ANSI/AAMI ST79:2006 and A1:2008, A2:2009 (Consolidated Text, Excerpt from Section 7.5.6).
3. "Comparative Efficacy of Hand Hygiene Agents", Sickbert-Bennett, Rutala, et al, Univ. North Carolina, Assoc. for Professionals in Infection Control (APIC), 2005.